



ASPEN IMAGING
HEALTHCARE

Europa60

Portable X-ray System Unit



GAME-CHANGING PORTABLE RADIOGRAPHY

Europa60 is a solution that adapts to your workspace, regardless of its size. With a compact design, its wireless and battery-operated system allows operators to capture X-ray images without requiring the dedicated physical space demanded by conventional radiographic systems.

Unlike other portable X-ray machines, Europa60 can acquire images directly on the device, allowing operators to review them first and then wirelessly transmit them to a PC via the built-in workstation.

DEVICE FEATURES

TWO USB PORTS

Used for powering peripheral devices and the laser component

7" TOUCH SCREEN DISPLAY

High Definition 7" touch screen display

ETHERNET PORT

Used for connecting exposure hand switch

CHARGING PORT

Simple battery charging - plugs directly into a residential 120 V outlet

HAND STRAP

Hand Strap for ease of use, supporting Europa's 5.5 lb lightweight



WHAT'S INCLUDED

- Europa60
- Carrying case
- Handswitch
- Charger
- Stylus pen
- Carrying strap



ALL-IN-ONE PC WORKSTATION

Integrated INTEL PC with
Windows 10+ compatible,
preloaded with DR Software



EXTREME LIGHT WEIGHT

5.5 lbs for ultimate
portability and convenience



REMOTE SUPPORT



Pre-loaded Remote Assistance
Software (AspenView)

SPECIFICATIONS

Europa60

POWER OUTPUT	120 W
kV / mA	60 kV(Fixed) / 2 mA(Fixed)
EXPOSURE TIME	0.03 ~ 1.3 sec
COLLIMATOR	Multi-leaf collimator
ANODE TYPE	Stationary anode
ANODE ANGLE	20°
FOCAL SPOT	0.8 mm
HEAT CAPACITY	6,000 J
BATTERY CAPACITY	22.2 VDC / 1,000 mAh
OPERATION	LCD touch screen / PC
WEIGHT	About 2.5 kg / 5.5 lb
DIMENSIONS	233 × 117 × 106 mm / 9.17" x 4.60" x 4.17"

HANDLED MODE

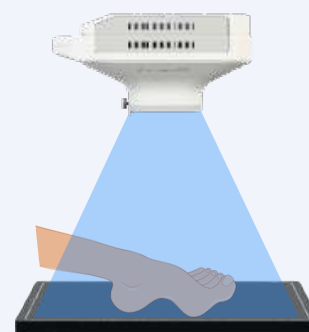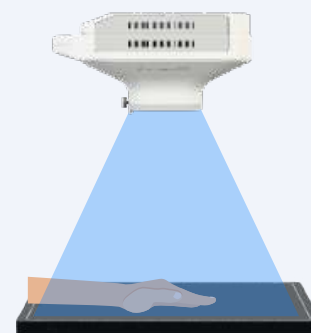



ASPEN IMAGING
HEALTHCARE

801 Jupiter Rd Suite 200, Plano, TX 75074, USA
sales@aspenimaging.com | (+1) 214-257-0113

WWW.ASPEN-IMAGING.COM

DISTRIBUTED BY

About Aspen Imaging

Aspen Imaging Healthcare is a Texas-based Original Equipment Manufacturer (OEM) providing complete imaging solutions for x-ray distributors and custom designed, private-labeled solutions for OEMs across the Americas.

Founded in 2015, Aspen Imaging has become one of the fastest growing Medical Imaging companies with offices in the United States, Mexico and South Korea.

